

ROS Packages for Mover4 and Mover6 Robot Arms

Contact: www.cpr-robots.com
crm@cpr-robots.com

1. Summary


The robot package (cpr_mover) and the RViz plugin (cpr_rviz_plugin) allow to integrate the Mover4 and Mover6 robot arms into ROS environments. An integrated JointTrajectoryAction server allows to perform MoveIt generated motions.

2 Tested Environments:

- ROS Indigo, Ubuntu 14.04 LTS


3 Architecture

The standard architecture consists of the core, RViz and the cpr_mover node.

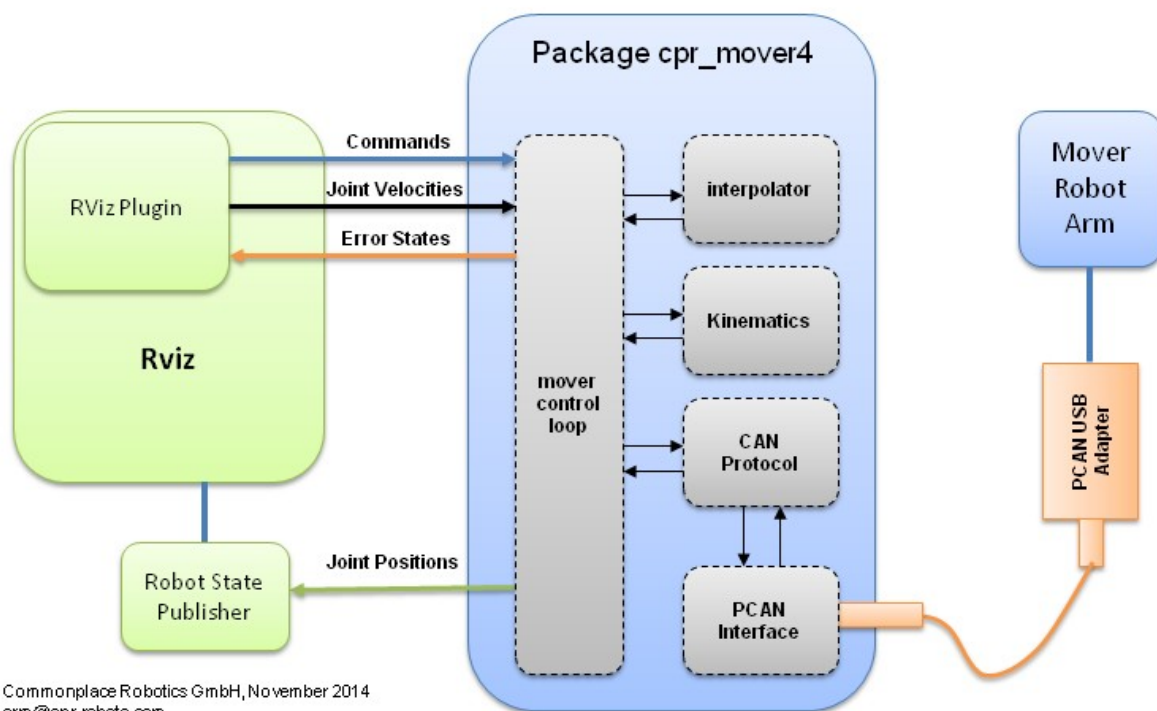

Commonplace Robotics GmbH, November 2014
crm@cpr-robots.com

The interaction with MoveIt is shown in the following slide.

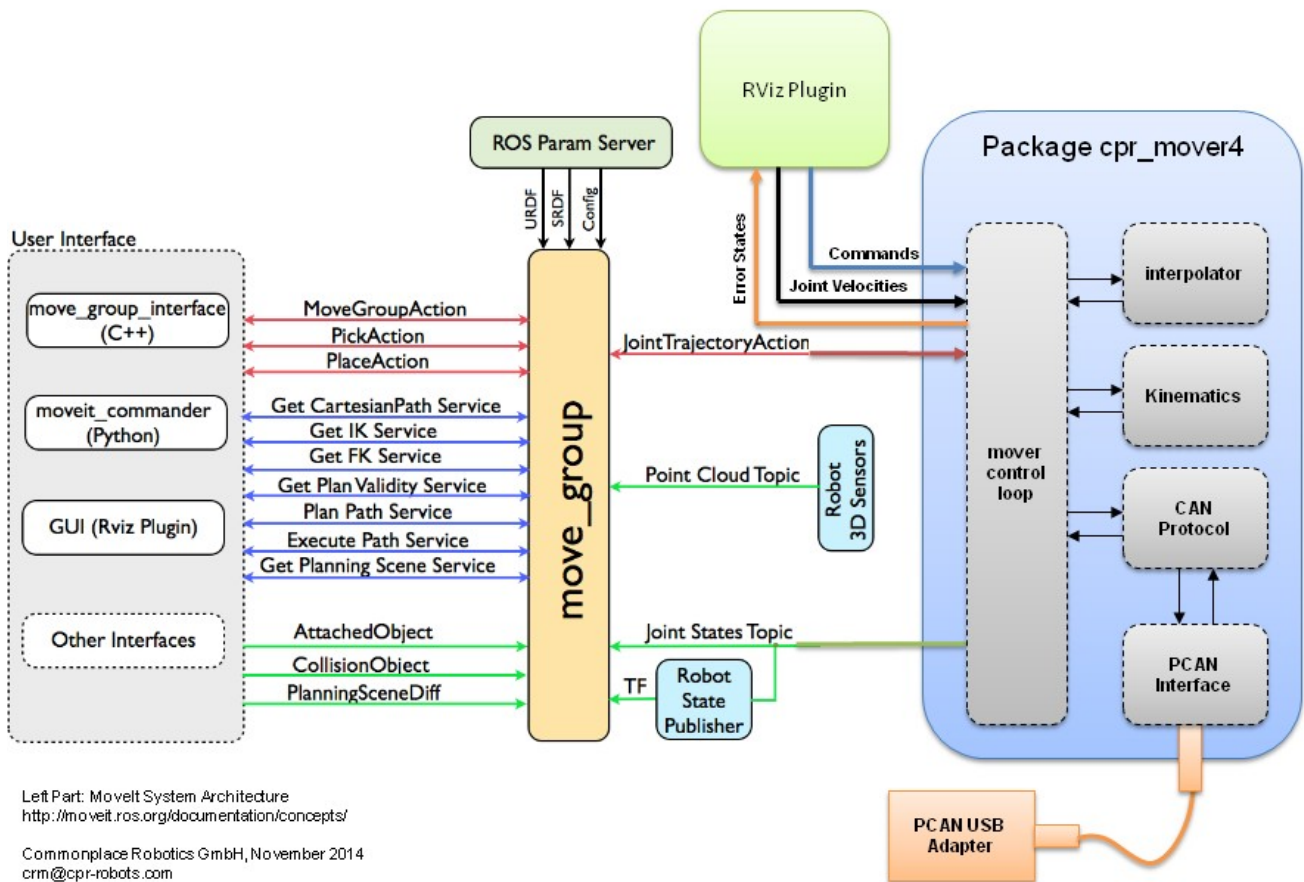

2. Download & Installation

2 Package cpr_mover

Download at http://www.github.com/CPR-Robots/cpr_mover

Save in your catkin-src directory, e.g. `~/catkin_ws/src/cpr_mover`

Compile with `catkin_make`

3 Plugin cpr_rviz_plugin

Download at http://www.github.com/CPR-Robots/cpr_rviz_plugin

Save in your ROS workspace, e.g. `~/ros-workspace/cpr_rviz_plugin`

Compile with `rosmake`

4 Peak PCAN USB-Adapter

The robot arm is connected by a Peak System PCAN-USB adapter. To use the adapter the according driver needs to be installed:

- Download the current driver package from Peak

- Extract driver package
- change into the directory, e.g. `cd peak-linux-driver-7.10`
- `make clean`
- We want to install the chardev version of the driver: `make NET=NO_NETDEV_SUPPORT`
- `sudo make install`
- `reboot`
- the command `cat /proc/pcan` should show:

```
*----- PEAK-System CAN interfaces (www.peak-system.com) -----
*----- Release_20140121_n (7.10.0) Mar  3 2014 14:13:14 -----
*----- [mod] [isa] [pci] [dng] [par] [usb] [pcc] -----
*----- 1 interfaces @ major 249 found -----
*n -type- ndev --base-- irq --btr- --read-- --write- --irqs-- -errors- status
32  usb -NA- ffffffff 255 0x001c 00000000 00000000 0000000c 00000000 0x0000
```

Especially below the `ndev` word there should be `-NA-`


When the adapter is recognized by the OS the red LED on the adapter is on. When the adapter is connected by the software, it blinks slowly. When data is transmitted it blinks faster.

3. Start

- Start `roscore`
- Open a terminal and `cd` into the plugin directory: `cd ~/ros_workspace/cpr_rviz_plugin/`
- Start Rviz with the according launch file: `roslaunch cpr_mover4.launch`
or `roslaunch cpr_mover6.launch`
- The CPR plugin should be loaded, otherwise load with the menu entry.
- Open a second terminal and `cd` into the `cpr_mover` directory: `cd ~/catkin_ws/src/cpr_mover`
- Start the node: `roslaunch cpr_mover cpr_mover`
It is important to first start Rviz, because the launch file sets the `robot_type` parameter to choose between Mover4 and Mover6, `cpr_mover` uses this parameter. If RViz is not used it can be set in another way to the parameter server.
- If the package does not start: `source ./devel/setup.bash`

Now move the simulated robot. If this on moves connect the real robot:

- Press “Connect” in the RViz plugin. The status changes from “not connected” to another value
- Press “Reset”: This button loads the hardware joint status, the 3D graphics adapts. Also the status changes
- Press “Enable”: Not the status changes to “0x00” or “No Error” and the robot can be moved now.

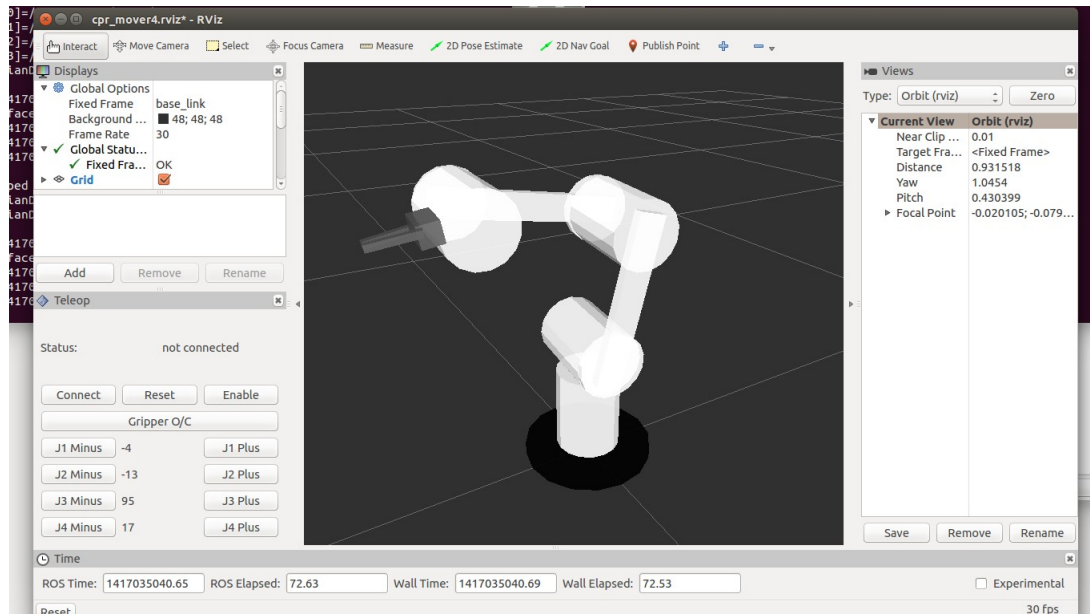

4. Interface Specification

Publishing Information takes place in `cpr_mover::CommunicationROS()` within the robot cycle time (standard: 20 Hz).

1 JointState Publisher

Type: `sensor_msgs::JointState`

Name: `/joint_states`

Provides 6 joint values + 2 gripper joint values of the robot in radian. The values are setpoint values, not the hardware values. Joint names are Joint0 to Joint5 and Gripper1, Gripper2.

The gripper is not commanded as servo joint, but using digital out by commands, see below.

2 Error Code Publisher

Type: `std_msgs::String`

Name: `/CPRMoverErrorCodes`

Provides a string with the current status of the robot arms hardware joints, the error codes.

3 Joint Velocity Subscriber

Type: `sensor_msgs::JointState`

Name: `/CPRMoverJointVel`

When there are no points to replay from the actionServer, the robot reacts to the jog values in these messages. The values in `msg->velocities[]` are percent values with respect to the `maxJointVelocity` defined in the source code. The allowed range is `[-100.0 .. 100.0]`.

4 Commands Subscriber

Type: `std_msgs::String`

Name: `/CPRMoverCommands`

The Mover reacts to commands received with these messages. Commands are:

- Connect
- Reset
- Enable
- GripperOpen
- GripperClose
- Override ppp where ppp is the percent value of the override (integer)

5 TrajectoryActionServer

Type: `actionlib::SimpleActionServer<control_msgs::FollowJointTrajectoryAction>`

Name: `cpr_mover/follow_joint_trajectory`

6 GripperCommandActionServer

Type: `actionlib::SimpleActionServer<control_msgs::GripperCommandAction>`

Name: `cpr_mover/gripper_command`

7 Trouble Shooting

- Slow, interrupted motions of the robot: Check if there are two TeleopPanels active in RViz. If this is the case delete both and load one again.